


Flat Specimen Holder

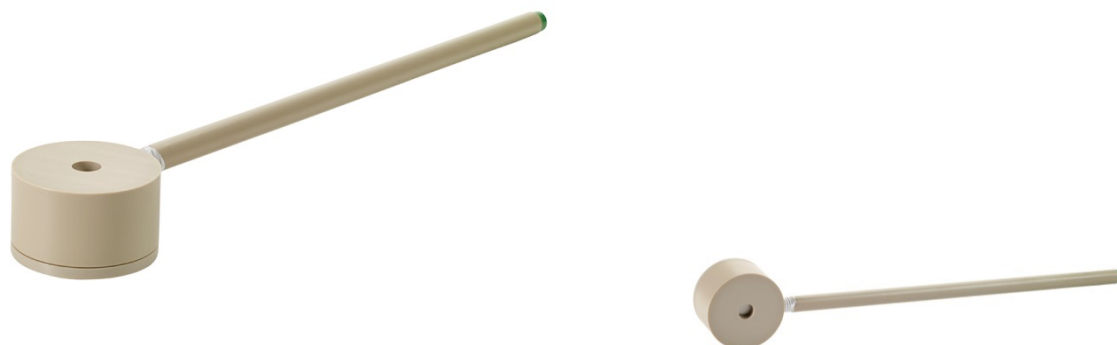

Holder Description

The Flat Specimen Holder (part number 990-00403) has been designed to hold flat, circular or square samples for use in the Multiport Corrosion Cell. This optional holder is designed to handle circular samples from 25 to 30 mm in diameter or square samples with sides length up to 23 mm. The actual exposed electrode face is 10 mm in diameter with an area of 0.785 cm². Samples can be up to 7 mm thick.

Made from tough, durable PEEK, the holder can withstand temperatures up to 80 C. Contact is made on the front of the sample through the use of Pogo pins.



Flat Specimen Holder in the Multiport Corrosion Cell.

*Patent Pending - WO 2016/016214 A1.